

The Rath Mechanical

TREATMENT TABLE

for the McKenzie technique



"Adjustable height with flexion and extension at the touch of a foot pedal."

Hill Laboratories Company, together with Wayne Rath P.T., developed the adjustable-height Rath Mechanical Treatment Table specifically for use with the well-known McKenzie Technique.

The Rath Mechanical makes flexion or extension positioning

of your patient easy with a touch of the electric foot pedal. Table movement can also be directly controlled by the patient with a convenient, 12-volt hand switch which neatly fastens under the table-top when not in use.

A rocker foot pedal allows the Rath Mechanical to be adjusted in height from 22" to 35" from either side of the table providing the ideal working height for

mobilization, therapy and examinations. Widths range from the standard 24" to 27" and 30". A durable telescoping base safely encloses all mechanical parts.

The Rath Mechanical Treatment Table has a one-year warranty on all parts and labor.


Extension Position



Flexion Position




~ See other side for featured options and specifications ~

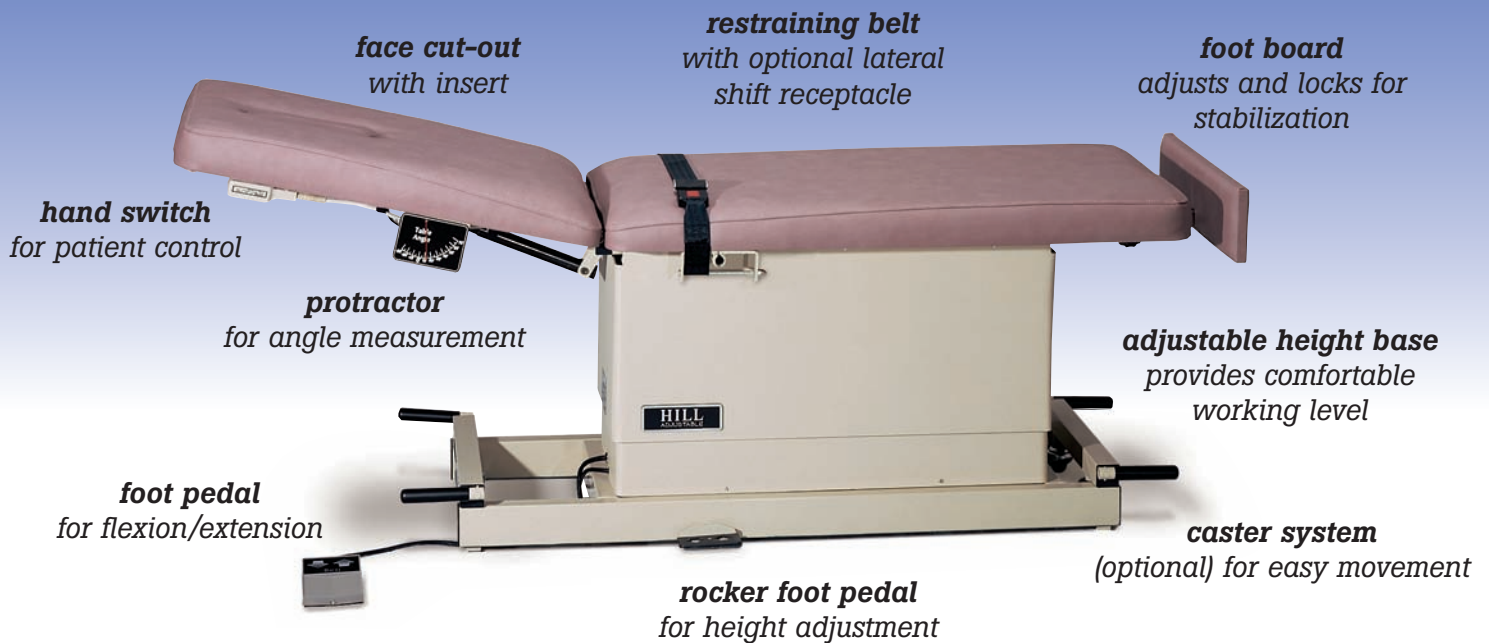

The Lateral Shift Device


SPECIAL OPTION

The lateral shift accessory for the Rath Mechanical Treatment Table enables the clinician to sustain a side gliding position of the lumbar spine while exploring

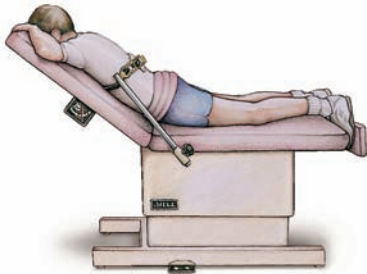

the effect of sustained and/or repeated passive extension procedures. The patient's response is monitored closely, looking for the presence of centralization or peripheralization of the signs and symptoms. The shifter can be of benefit in the management of the patient with a relevant lateral shift or lateral compartment who is unable to achieve full benefit with self-generated forces.

These procedures should only be used after the clinician has received proper training in Mechanical Diagnosis and Therapy of the lumbar spine.

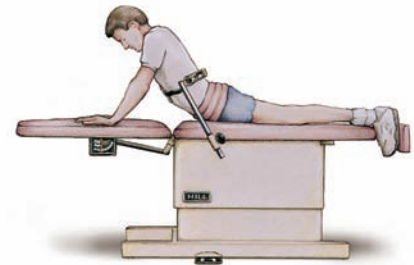

Specifications

- Elevates from low to high in 16 seconds
- Thermo-plastic base in beige or gray
- Table length-76"
- 2" thick Ultra-Cell Foam™ for comfort and shape retention
- Durable vinyl upholstery in choice of 20 colors
- Lifting capability - 500 pounds
- Electric requirement - 110 volts, 60 Hz, 4.5 amps
- 12 volt hand switch
- One-year warranty on parts and labor
- Shipping weight approximately 225 pounds

ForscoLabs
Exclusive Distributor - HILL Laboratories